


Configuring mSupply Mobile on Device


Initialising the store

After installing the mSupply Mobile app, you must do an initialisation of the data for the store that the tablet is going to be used in. On opening the mSupply Mobile app, the first screen looks as follows:

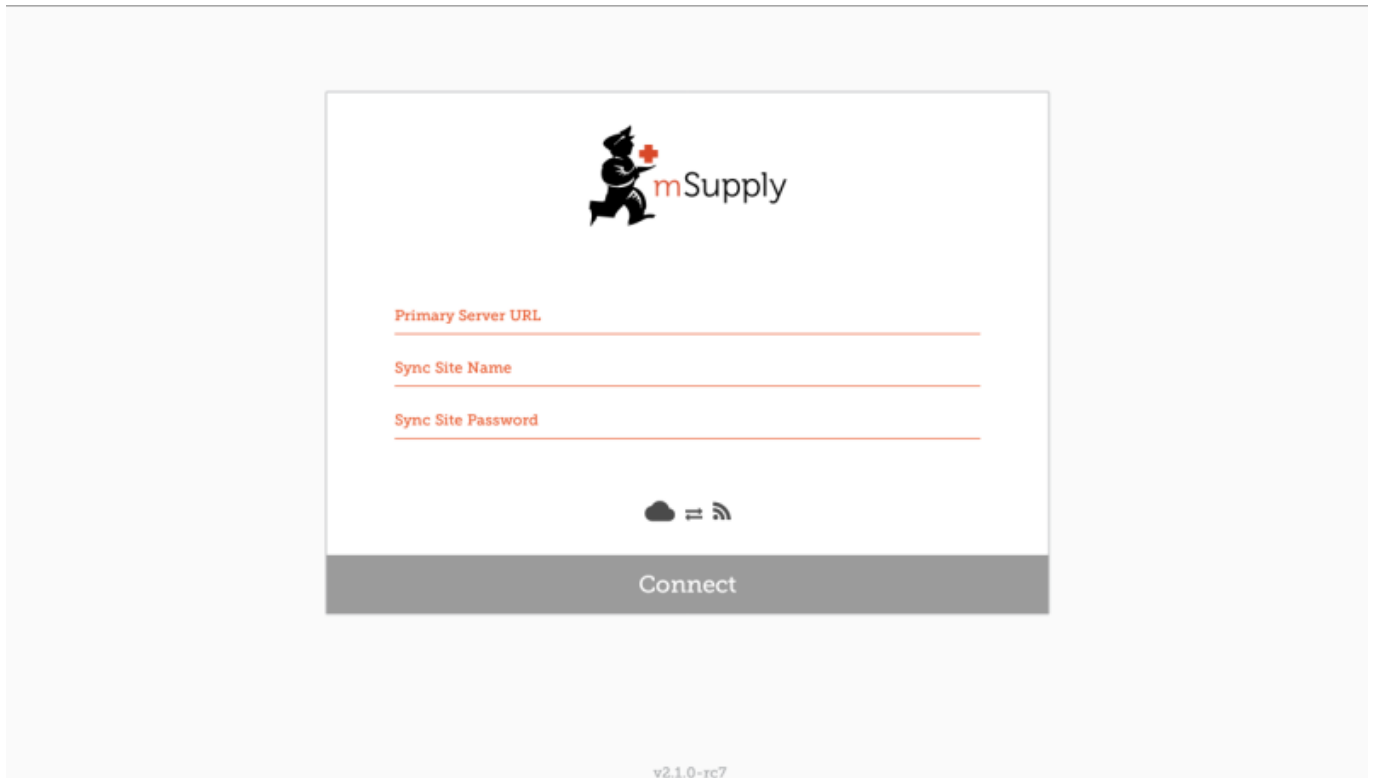

You must enter:


- The correct URL for your mSupply server, usually using SSL/HTTPS protocol for security, e.g <https://yoursite.msupply.org>
- Sync site name as configured on the Primary Server.
- Password for the sync site name as configured on the Primary Server



https://my.server.org:(mySSLportNumber)


My clinic's sync site name

.....

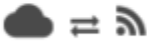

Connect


Once you have entered the information, click "Connect" to start sync:



http://192.168.4.16:8088

demo_clinic


....



Initialising...

Initialisation can take between 1 and 15 minutes, depending on the number of records on the server related to the site being synced. Once the server has found all the records to sync, mSupply Mobile will start pulling them adaptively between 20 and 500 records at a time depending on the connection speed. If connection is lost or you closed the app, you can resume sync where it stopped by entering

the exact same details. It will start again from scratch if you get it wrong before getting it right again!

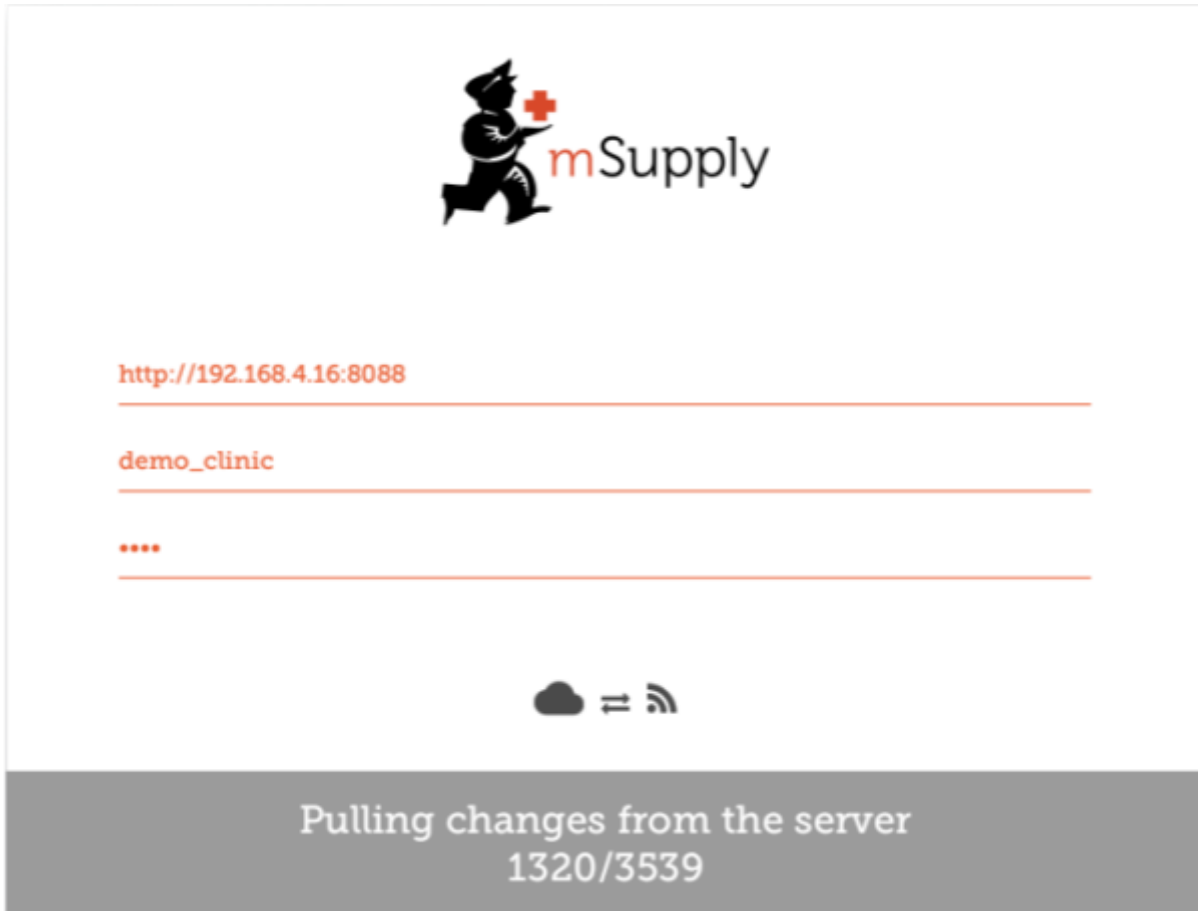


Logging in to the store

Once complete, you will be taken to the login screen. For guidance beyond here see the [mSupply Mobile usage guide](#).

Limiting Initial Sync Record Count

If the amount of records to sync is large (say around 50,000) it may take hours to sync, particularly on a *poor* internet connection. For this reason it is possible in mSupply to limit the amount of records up to a specified amount of months ago.

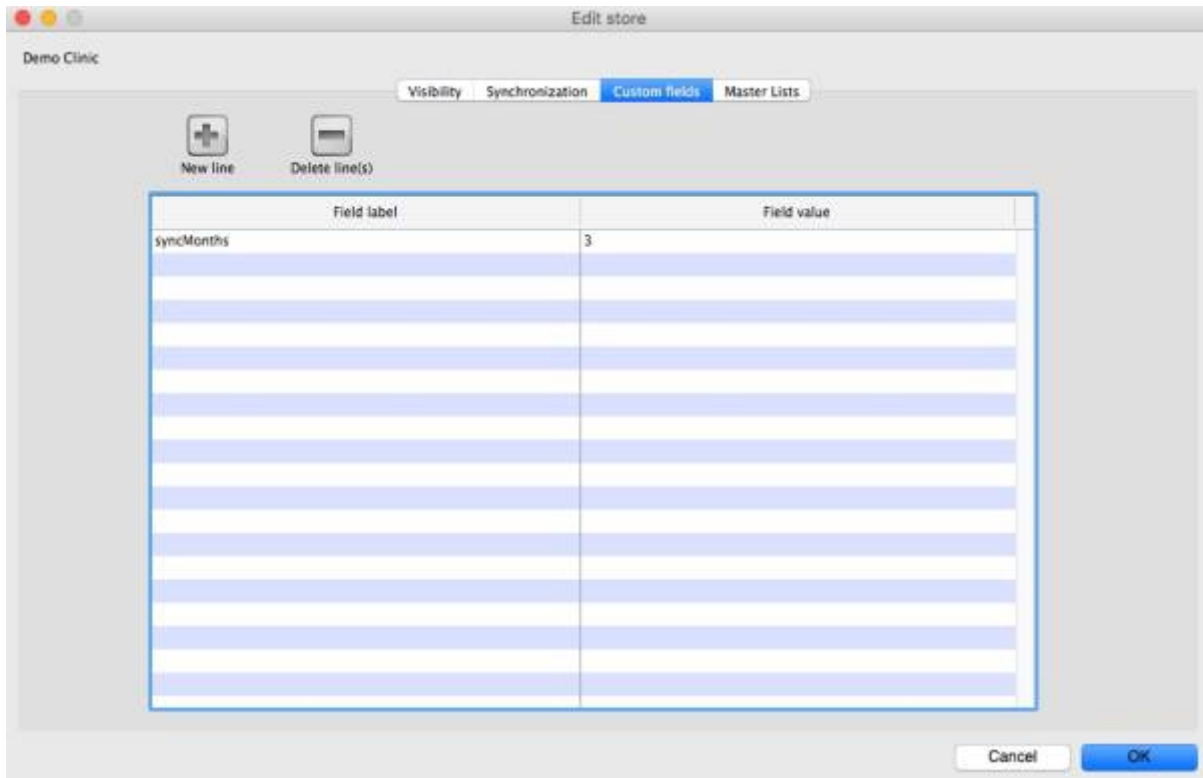

1. You need to change the store synchronisation type from “collector” to “active” or “none” to unlock adding custom data
2. Add a custom data row with “Field label” `syncMonths` and the numeric value desired in “Field value”
3. Change the store synchronisation type back to collector

Now when you do initial sync, the store will only get transactions and stock takes for as far back as specified from the current date.



Related Links

- [mSupply mobile user guide](#)

Return to: [mSupply Mobile Home page](#)

From: <https://wiki.msupply.foundation/> - **mSupply Foundation Documentation**

Permanent link: https://wiki.msupply.foundation/en:mobile:setup_device:1?rev=1575351436

Last update: **2019/12/03 05:37**

